


PFA High Purity Sample Introduction System for Agilent ICP-MS systems

Inert sample introduction technology for improved background and detection limits

An easy to install, inert, high purity, o-ring free upgrade kit which includes two PFA MicroFlow nebulizers, PureCap endcap, 35mm PureChamber PFA spray chamber, adapter, o-ring free injector and precision quartz torch. A reliable ultra pure analysis technology to complement Agilent ICP-MS systems

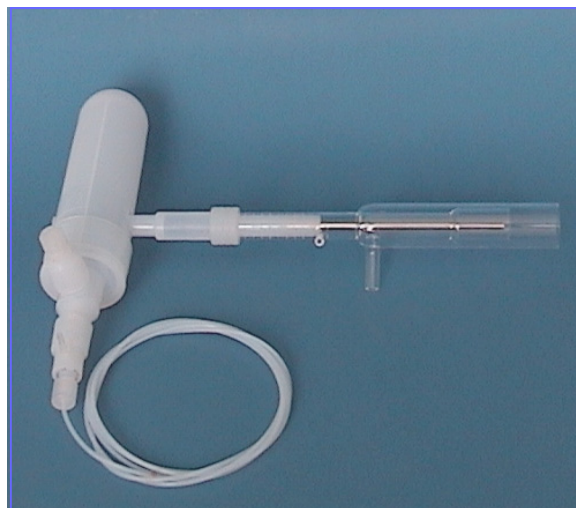

All Teflon, o-ring-free full upgrade kit for the Agilent ICP-MS system.

- Easy installation on all Agilent ICP-MS systems.
- Clean Teflon® PFA components for lower background equivalent concentration (BEC) and detection limits.
- O-ring free connections —ensures lowest BEC's, highest chemical resistance and easy cleaning.
- Chemically resistance—suitable for nearly all samples, including strong acids, alkalis and organic solvents.
- Sapphire or platinum injector.
- Two self-aspirated or pumped PFA nebulizers included.
- Low sample consumption—ideal for VPD and preconcentrated samples

